

FRED® EASYPORT® PLUS


Small and powerful defibrillator for anyone - anytime - anywhere



FRED® EASYPORT® PLUS

The ideal companion on the go

FRED® easyport® plus, the pocket-size defibrillator, fits easily into every doctor's bag and emergency kit. It is the ideal companion during any intervention.

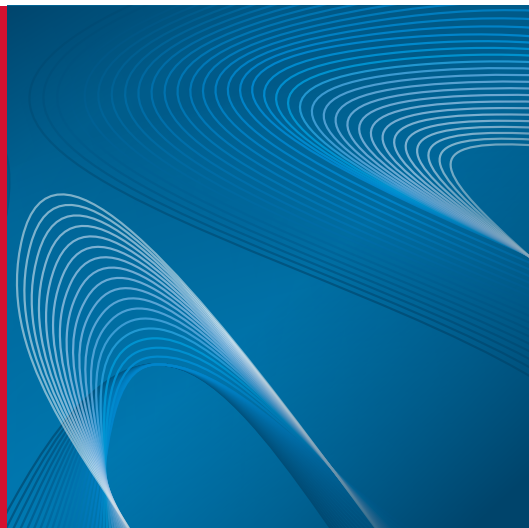

IT CONCERNS ALL OF US

Every year, more than 7 million people fall victim to a sudden cardiac arrest worldwide. A sudden cardiac arrest can hit anyone, anywhere and at any time.

Defibrillation delivers electrical energy to the heart, restoring a proper heart rhythm. A defibrillator in combination with well-performed chest compressions is the only therapy that can save a patient's life.

EASY TO CARRY

As a pocket-size defibrillator, the FRED easyport plus weighs only 780 grams and combines small size with great functionalities.

SHOCK AS SOON AS IT IS REQUIRED

International guidelines emphasise the importance of minimising the delay in shock delivery.

To reduce this delay, the FRED easyport plus has a quick analysis system and loads the required energy systematically so that the shock can be delivered as soon as the device recommends it.

ADAPTED TO ANY USER

The FRED easyport plus is suitable for frequent as well as non frequent users. It can be operated with a rechargeable or a disposable battery. The external charger can charge two batteries simultaneously.

CPR FEEDBACK

Defibrillation alone is not enough. To maintain a good supply of oxygen to the organs, cardiopulmonary resuscitation (CPR) needs to be performed correctly to restore a proper blood flow.

The FRED easyport plus helps you to achieve highest efficiency in chest compressions by:

- ❖ indicating the correct depth, frequency and recoil*
- ❖ setting a regular rhythm through the metronome

When using the ARGUS LifePoint sensor, the rescuer benefits from audio and visual assistance through on-screen colour-coded messages while performing chest compressions.

*Recoil not available in France, Germany, UK and USA



INFORMATION TRACKING

The FRED easyport plus records intervention data, allowing doctors and emergency physicians to get a better overview of the patient's situation.

- ❖ CPR quality
- ❖ ECG curve
- ❖ Impedance curve
- ❖ Events

INFORMATION RETRIEVAL

Thanks to USB and Bluetooth functionalities, saved files can easily be reviewed and managed by SCHILLER's data management solution SEMA.

FIND THE DEVICE THAT SUITS YOU

The FRED easyport plus is available in three versions to fit the different user levels.



AED FIRST

The FRED easyport plus semi-automatic AED with its grey cover requires the rescuer's confirmation to deliver a shock.

An ECG curve can be displayed on the large LCD colour screen.

Suitable for first responders.



AED AUTO

The fully automatic FRED easyport plus does not require any action from the rescuer. After the heart rate analysis, the device determines whether a shock is necessary. If it is, the shock is delivered automatically.

An ECG curve can be displayed on the large LCD colour screen.

Suitable for trained lay responders.



AED MANUAL

Distinguishable by its red cover, the manual FRED easyport plus allows doctors and paramedics to switch to manual override. Thanks to the ECG display, healthcare professionals can determine if a shock is necessary.

Suitable for physicians.



SCHILLER
FRED easyport plus

00:00:04
PLUG AND APPLY ELECTRODES



ELECTRODES
FRED PA-1,
Adult, 1000mAh Pro Life Case 1,
Type 2, CARDIOMIT AT-101, 1000mAh



SMALL AND POWERFUL

The FRED easyport plus weighs only 780 grams and delivers shocks of up to 200 joules. The Multipulse Biowave® technology ensures an effective and safe defibrillation shock while minimising damage to the myocardium.

GAIN PRECIOUS SECONDS

The defibrillation electrodes are pre-connected for a quick application on the patient's chest. Heart rate analysis starts right away to reduce time until the shock is delivered.

ALWAYS READY

The FRED easyport plus performs regular automatic self-tests to ensure it works correctly. Its status is indicated by the flashing on/off button.

DATA TRANSMISSION TO SEMA

Intervention data can be sent via USB or Bluetooth for review by healthcare professionals.

DESIGNED TO LAST

With IP44 classification, the FRED easyport plus resists entry of splashing water.





POCKET-SIZE DEFIBRILLATOR

The FRED easyport plus measures only 150 x 143 millimeters – this photo shows its original size (1:1)!

SUITABLE FOR ALL ENVIRONMENTS

The FRED easyport plus has a large colour LCD screen which automatically adjusts its brightness to help visualise indications at a glance.

INTUITIVE PAEDIATRIC MODE

While using the same electrodes, just press this button to change the patient mode.

CPR FEEDBACK

The integrated metronome sets the correct rhythm and guides the rescuer during the intervention.

In addition, the ARGUS LifePoint sensor assists the rescuer with audible and visual cues on screen to ensure effective chest compressions.

OPTIMISED USER INTERFACE

Three simple steps with illustrations and spoken instructions for a clear understanding.



Americas
SCHILLER Americas Inc.
Doral, Florida 33172

North America:
Phone +1 786 845 06 20
Fax +1 786 845 06 02
sales@schilleramericas.com
www.schilleramericas.com

Latin America & Caribbean
Phone +1 305 591 11 21
Fax +1 786 845 06 02
sales@schilleramericas.com
www.schilleramericas.com



Asia
SCHILLER Asia-Pacific / Malaysia
52200 Kuala Lumpur, Malaysia
Phone +603 6272 3033
sales@schiller.com.my
www.schiller-asia.com



Austria
SCHILLER Handelsgesellschaft m.b.H.
A-4040 Linz
Phone +43 732 70 99 0
Fax +43 732 757 000
sales@schiller.at
www.schiller.at



China
Alfred Schiller (Beijing) Medical Equipment
Co.,Ltd.
100015 Beijing, China
Phone +86 010 52007020
Fax +86 010 52007020-8016
info@schillerchina.com
www.schillermedical.cn



Croatia
SCHILLER d.o.o.
10000 Zagreb
Phone +385 1 309 66 59
Fax +385 1 309 66 60
info@schillerzg.hr
www.schiller.hr



France
SCHILLER Médical
F-67160 Wissembourg
Phone +33 3 88 63 36 00
Fax +33 3 88 63 36 49
info@schiller.fr
www.schiller-medical.com



France (distribution France)
SCHILLER France S.A.S.
F-77608 Bussy St Georges
Phone +33 1 64 66 50 00
Fax +33 1 64 66 50 10
contact@schillerfrance.fr
www.schiller-france.com



Germany
SCHILLER Medizintechnik GmbH
D-85622 Feldkirchen b. München
Phone +49 89 62 99 81 0
Fax +49 89 609 50 90
info@schillermed.de
www.schillermed.de



India
SCHILLER Healthcare India Pvt. Ltd.
Mumbai – 400 059, India
Phone +91 22 6152 3333 / 2920 9141
Fax +91 22 2920 9142
sales@schillerindia.com
www.schillerindia.com



Poland
SCHILLER Poland Sp. z o.o.
PL-02-729 Warszawa
Phone +48 22 843 20 89 / +48 22 647 35 90
Fax +48 22 843 20 89
schiller@schiller.pl
www.schiller.pl



Russia & C.I.S.
AO SCHILLER.RU
125040 Moscow, Russia
Phone +7 (495) 970 11 33
Fax +7 (495) 956 29 10
mail@schiller.ru
www.schiller.ru



Serbia
SCHILLER d.o.o.
11010 Beograd
Phone +381 11 39 79 508
Fax +381 11 39 79 518
info@schiller.rs
www.schiller.rs



Slovenia
SCHILLER d.o.o.
2310 Slovenska Bistrica
Phone +386 2 843 00 56
Fax +386 2 843 00 57
info@schiller.si
www.schiller.si



Spain
SCHILLER ESPAÑA, S.A.
ES-28232 Las Rozas/Madrid
Phone +34 91 713 01 76
Fax +34 91 355 79 33
schiller@schiller.es
www.schiller.es



Switzerland
SCHILLER-Reomed AG
CH-8912 Obfelden
Phone +41 44 744 30 00
Fax +41 44 740 37 10
sales@schiller-reomed.ch
www.schiller-reomed.ch



Turkey
SCHILLER TÜRKIYE
Okmeydani-Sisli – Istanbul
Phone +90 212 210 8681 (pbx)
Fax +90 212 210 8684
info@schiller.com.tr
www.schiller-turkiye.com



United Kingdom
SCHILLER UK Ltd.
Coatbridge, ML5 2DA
Phone +44 1619 072 555
sales@schilleruk.com
www.schilleruk.com

Device availability in your market is subject to regulatory approval. Subject to changes.

2021/01, EN



Article no.: 2.500782



SCHILLER
The Art of Saving Lives

Headquarters: SCHILLER AG, Altgasse 68, CH-6341 Baar, Switzerland, Phone +41 41 766 42 42, Fax +41 41 761 08 80, sales@schiller.ch, www.schiller.ch